


**Zirkonzahn®**

*Human Zirconium Technology*



## **THE ART OF COLOURING ZIRCONIA**

*Colouring liquids for your Prettau® Zirconia masterpieces*

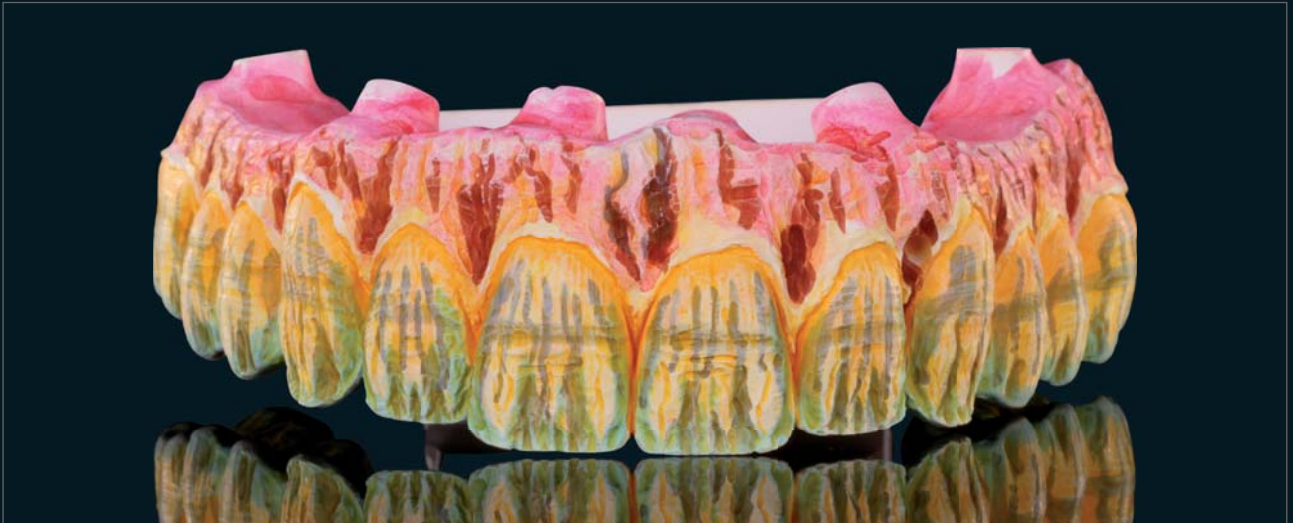

# OCCLUSALLY SCREW-RETAINED PRETTAU® BRIDGES INFILTRATED WITH COLOURING LIQUIDS USING ZIRKONZAHN'S SPECIAL TECHNIQUE

1. *Acquisition of the patient's situation with PlaneSystem® technology*
2. *Milling the titanium bar and the Prettau® Bridge with the compact M1 Wet Heavy Metal milling unit*
3. *Infiltrating with Zirkozahn's acid-free, water-based colouring liquids: with 100% burn-out biocompatible pigments, the liquids assist in realistic shade distribution and colour grading allowing for an accurate and individual characterisation*
4. *Sintering in the Zirkonofen 700: the bio-pigments are cleanly burnt out and the restoration appears in the desired tooth shade*
5. *Layering with ceramics, characterisation with ICE Zirkon 3D Stains by Enrico Steger and glazing*
6. *Anodising the titanium bar and bases to a golden colour to reduce the restoration's grey value*
7. *Bonding of the titanium bar and basis to the zirconia structures*
8. *The masterpiece is ready!*



**Zirkonzahn®**



**UPGRADE  
YOUR  
KNOWLEDGE!**

**JOIN OUR LECTURES  
IN THE UK FOR DENTAL  
TECHNICIANS AND  
DENTISTS!**



*Scan the QR  
code to see  
all upcoming  
conferences*

*Have a look at the complete case  
gallery at [www.zirkonzahn.com](http://www.zirkonzahn.com)*

## **HUMAN ZIRCONIUM TECHNOLOGY**

*Zirkonzahn Worldwide – Tel +39 0474 066 680 – [info@zirkonzahn.com](mailto:info@zirkonzahn.com) – [www.zirkonzahn.com](http://www.zirkonzahn.com)*



**MILLED WITH THE  
M1 WET HEAVY METAL  
CAD/CAM SYSTEM**



WEAA5498=

*All information is subject to change. Errors and omissions excepted. Version: 08/12/2017*